WRINKLE REDUCER

N Ū F A C E[®]

USER MANUAL

ENSURE YOUR BEST RESULTS!

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NUFACE® WRINKLE REDUCER

TRINITY+ ATTACHMENT

Note: Refer to NuFACE Trinity+ User Manual for more device Information.

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2) Magnetic – NuFACE Trinity+ Wrinkle Reducer Attachment

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The Trinity+ Wrinkle Reducer is a hand-held phototherapy device that emits Amber, Red and Infrared (IR) light energy via a light attachment head. The Amber, Red, and IR light energy is used in the treatment of full-face fine lines and wrinkles, and increase in local circulation within the perioral region (the area around the mouth).

IMPORTANT INFORMATION

- Use the Trinity+ Wrinkle Reducer attachment as described within this User Manual.
- The NuFACE Trinity+ device User Manual may also be referenced.
- The risks and dangers of using the attachment in any way that does not follow this User Manual are unknown.
- The Trinity+ Wrinkle Reducer was not tested for use over the eye socket or eyelid so the risks are unknown.
- Your Wrinkle Reducer attachment has a unique serial number located on the underside of the attachment.
- Register your new device at www.mynuface.com/pages/register.

CONTRAINDICATIONS FOR USE

- **DO NOT** use if you are pregnant, subject to epilepsy, or have photosensitive epilepsy.
- **DO NOT** use on the neck if diagnosed with or suspect a thyroid problem.
- **DO NOT** use if your skin is sensitive to light. People can be sensitive to light and not be aware of it. Certain drugs or cosmetics can cause skin sensitivity.
- **DO NOT** use over swollen, infected, or inflamed areas or skin eruptions.
- **DO NOT** treat the skin more than once in a 24 hour period. Overuse may result in dryness of the skin.
- **DO NOT** use over topical or systemic steroids.
- **DO NOT** use on infants or children under the age of 14.
- * Consult your healthcare provider if any of the contraindications above apply to you.

WARNINGS

- If you are under the care of a physician, consult your physician before using this device.
- This device is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Apply stimulation only to normal, intact, clean healthy skin. **DO NOT** use over swollen, infected, or inflamed areas or skin eruptions (e.g., phlebitis, thrombophlebitis, broken capillaries, varicose veins, etc).

Consult your healthcare provider:

- If you have any pre-existing or potential health conditions that may be impacted by use of the device.
- If you have undiagnosed or known pain syndromes or painful areas.
- If you have any underlying skin conditions.
- If you have undiagnosed or known pain syndromes.
- **DO NOT** operate Trinity+ Wrinkle Reducer in the presence of electronic monitoring equipment (e.g., heart monitors, heart tracing or in close proximity (e.g. 1 m) to shortwave, microwave or electromagnetic therapy equipment as it may cause the medical equipment to malfunction. Such equipment may produce instability in the Trinity+ Wrinkle Reducer.

WARNINGS (CONTINUED)

- If you suspect the device is overheating, turn the device OFF immediately.
- **DO NOT** use flammable aerosol products or oxygen near the device.
- **DO NOT** look directly into the Trinity+ Wrinkle Reducer LEDs or shine in the eyes as this may cause harm.
- Overuse of this device may result in dryness of the skin.
- Certain drugs or cosmetics can cause skin sensitivity.
- If device is not functioning properly due to damage caused by dropping or submerging in water.
- If you suspect the device to be malfunctioning, or not operating properly please contact NuFACE.
- Never drop or insert any object into any openings on the device.
- Device should not be used at extreme temperatures.
- Keep all electrical appliances (including the Trinity+ Wrinkle Reducer) away from water (including baths, showers, and sinks).
- **DO NOT** use while driving, operating machinery, or during any activity which may put the user or the public at unnecessary risk of injury.
- Use of controls or adjustments, or the performance of procedures other than those specified herein, may result in injury.
- **DO NOT** disassemble the attachment as this may cause damage, malfunction, electrical shock, fire or injury. There are no user-serviceable parts inside.
- **DO NOT** use over swollen, infected, or inflamed areas or skin eruptions. Do not treat the skin more than once in a 24 hour period -overuse may result in dryness of the skin. If you have any skin condition concerns, please consult your healthcare provider before use.

Discontinue use of the device and contact your healthcare provider if you experience any of the following:

- Persistent pain or discomfort.
- Negative or harmful reaction(s) or injury.
- Symptoms of hypertension (dizziness or light-headedness).
- Skin irritation or burning.
- Unpleasant sensations on skin during use.

Note: Disassembly of the NuFACE Trinity+ devices and/or attachments will void the warranty.

PRECAUTIONS

- The Trinity+ Wrinkle Reducer is designed for cosmetic use only and for individuals in good health. If this is not the case for you, please consult your healthcare provider before use.
- If using the Trinity+ Wrinkle Reducer causes persistent pain or discomfort, discontinue use immediately.
- The Trinity+ Wrinkle Reducer is intended to be stored and used at room temperature. Do not store Trinity+ Wrinkle Reducer in the sun or on a hot surface. High temperatures may be damaging. The device should be operated, transported and stored in the manner specified in this manual. Failure to use and maintain the device in accordance with the instructions outlined in this manual will void the warranty.

PRECAUTIONS (CONTINUED)

The NuFACE Trinity+ Wrinkle Reducer is intended for cosmetic use on the face as described. Any consequence resulting from misuse, application to other parts of the body or face, connection to improper voltage sources, dirty surface, or any other improper applications is neither the responsibility of NuFACE nor its affiliates.

BEFORE USING YOUR NEW ATTACHMENT

Follow these important steps before using your new attachment.

Follow steps 1-5 of the NuFACE Trinity+ Device (sold separately) Manual "BEFORE USING YOUR NEW DEVICE" section.

Note: The Trinity+ Wrinkle Reducer's unique technology utilizes flashing LEDs to address facial wrinkles. Begin treatment in a well-lit room if you find flashing lights bothersome.

1) READ the NuFACE Trinity+ Device User Manual (sold separately) Reference the section in the Trinity+ Device User Manual titled "BEFORE USING YOUR NEW DEVICE" before changing attachments.

2) CONNECT YOUR TRINITY+ WRINKLE REDUCER

The NuFACE Trinity+ device is designed with magnetic attachments.

- Remove the Facial Trainer attachment.
- Connect the Wrinkle Reducer attachment onto the NuFACE Trinity+ device.

The Trinity+ device will NOT turn on if the attachment is not securely in place.

3) TEST YOUR SKIN SENSITIVITY

IMPORTANT: Test your skin to make sure you are not sensitive to light.


- Place the Trinity+ Wrinkle Reducer on your forearm for **3 minutes** and then check the area.
- If your skin turns red and the redness lasts for more than 2 hours, you may be sensitive to light therapy and should consult with your healthcare provider before continuing to use the Trinity+ Wrinkle Reducer.

TREATMENT AREAS

The Trinity+ Wrinkle Reducer is designed to address large surface areas of the face. The following instructions are designed to provide users with the most effective results when followed correctly.

QUICK START

- 1) **Prep** – Start with clean, dry skin
- 2) **Treat** – Place the NuFACE Trinity+ with Wrinkle Reducer Attachment directly on skin of one of the 7 targeted areas. Hold in place until you hear the beep (3 minutes), then move to the next area. Do not treat upper eyelid.
- 3) **Finish** – Proceed to your NuFACE® microcurrent treatment or apply your favorite skincare products.

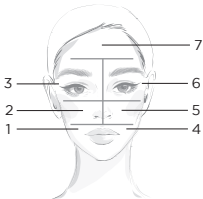HOW TO USE YOUR TRINITY+ WRINKLE REDUCER ATTACHMENT

Get ready to give wrinkles the red light.

PREPARE YOUR DEVICE AND ATTACHMENT

- Remove your device from the charging cradle.
- Ensure your Wrinkle Reducer attachment is securely in place and turn it ON by pressing the ON/OFF button.
- An ascending audible beep is emitted and the device's power lights will illuminate indicating that your device is ON.

DEVICE AUDIBLE SOUND	DEVICE INDICATION
Ascending Beeps & Hums	Powered ON
Descending Beeps	Powered OFF
Dual-Tone Beep	End of Glide Interval - 3 minutes (Tone occurs only during treatment)
Low-Tone Beep	Battery Needs To Be Charged (Tone occurs only during treatment)

HOW TO USE YOUR TRINITY+ WRINKLE REDUCER ATTACHMENT (CONTINUED)

1) PREP

- Cleanse your face and neck with a cleanser.

2) BEGIN TREATMENT

- Place the Trinity+ Wrinkle Reducer in direct contact with your skin, over the area to be treated.
- Your Trinity+ Wrinkle Reducer has a skin sensor that will illuminate the red LEDs when in complete contact with your skin.
- Press firmly enough to keep the Trinity+ Wrinkle Reducer in direct contact with your skin.
- The device should be used on lines and wrinkles within the section, not specifically where the number appears on the image.
- Treat each section for 3 minutes until the beep sounds.

3) FINISH TREATMENT

- Turn OFF the Trinity+ device after treatment.
- A descending audible beep is emitted notifying you that your device is turning OFF.
- The device will automatically turn off after 21 minutes.

If the Trinity+ Wrinkle Reducer is not touching the skin, your device will emit 3 beeps notifying you to hold directly on your skin.

TREATMENT FREQUENCY

- Use no more than 3 minutes per area, up to 21 minutes a day for the entire face. Do not overlap areas.
- For the first 60 days, use 5 days per week.
- After 60 days, use 2-3 days per week for maintenance.
- Allow 24 hours between treatments.
- For more advanced treatment options and holds, download the NuFACE App, or visit our website at www.NuFACE.com.

RESULTS

Results may vary based on skin type and skin condition. Best results are achieved with regular use.

NUFACE SMART APP

Download the NuFACE SMART APP for additional how-to instructions, and customized treatment options.

- You can download the NuFACE Smart App from the App store or the Google Play Store.
- To pair your device, in the App, click on **'Add Device'** and follow the on-screen instructions to connect. The NuFACE app is available in many countries around the world.

CARING FOR YOUR WRINKLE REDUCER ATTACHMENT

Cleaning

- Slightly dampen a soft, lint-free cloth with water to clean your device and attachment.
- **DO NOT** immerse device or treatment attachment in water at any time.

CARING FOR YOUR WRINKLE REDUCER ATTACHMENT (CONTINUED)

Storage

- The Trinity+ device and accessories are intended to be stored and used at room temperature. Use caution prior to using the device if it has been stored or not used for an extended period of time.

Cleaning

- Slightly dampen a soft, lint-free cloth with water or rubbing alcohol to clean your device and attachment. **Do not immerse device or treatment attachment in water at any time.**

Disposal

IMPORTANT: If disposing the accessories, follow local regulations for proper disposal.

TROUBLESHOOTING

Device Not Turning On

1) Check to see if your Trinity+ device is fully charged.

- If the device does not turn ON, the battery may be drained.
- A full charge cycle of ten (10) hours may be required for the device.
- Reference the 'NuFACE Trinity+ Device User Manual' for additional troubleshooting.

2) Check to see if the attachment is completely secured to the Trinity+ device.

IMPORTANT: The Trinity+ device will NOT turn on if the attachment is not securely in place.

DO NOT attempt to troubleshoot by disassembling the product as this will void the warranty and may cause injury.

For additional troubleshooting or information specific to the device, reference the NuFACE® Trinity+ Device User Manual.

CLIENT SUPPORT

If you are still experiencing issues with your device, contact NuFACE Client Services at **info@nuface.com** or call **888-360-9171** (US Only) between 8am PST to 5pm PST.

WARRANTY

NuFACE® warrants that all NuFACE® Trinity+ devices and attachments are free from manufacturer defects in materials and workmanship for a period of two (2) years from the date of your original purchase. Our obligation under this Warranty is limited to replacement of the product which will not extend the period of your warranty. Disassembly of the NuFACE Trinity+ devices and/or attachments will void the warranty. This Warranty does not cover damage covered by abuse, improper care, application of harmful chemicals, product that may have been altered in any way, loss of the device and product purchased from an unauthorized seller.

NO OTHER WARRANTY EXPRESSED OR IMPLIED, WHETHER OF FITNESS FOR A PARTICULAR PURPOSE OR OF MERCHANTABILITY OR OF ANY OTHER KIND, SHALL EXIST IN RESPECT TO SUCH PRODUCT.

WARRANTY (CONTINUED)

For Devices Purchased in the US: For more information about your warranty or to initiate a warranty claim and receive a return authorization number, contact our Client Services Team via our website at **www.NuFACE.com** or by calling **888-360-9171**. You must show a receipt from an authorized seller as proof of purchase.

For Devices Purchased Outside of the US, International Returns, and International Warranties: For all product purchases and returns outside of the United States, warranty and returns are per your place of purchase: spa, retailer, or distributor. Please consult your receipt or place of purchase for their return and warranty policy.

Conforms to EN IEC/UL: Device - 60601-1, 60601-1-6, 62366-1, 60601-1-11, 60601-2-57, 60529, 62471.

SYMBOLS

One or more of the following markings may appear on the Trinity+ Device, User Manual or packaging.

European Conformity

Comply with European Restriction of Hazardous Substances Directive 2002/95/EC

European Representative

Product must be disposed of in accordance with the WEEE directive

Part Number

FCC-approved equipment authorization mark

Refer to Instruction Manual

Class II double insulated electrical device

Manufacturer

Regulatory Compliance Mark

UK Conformity Assessed

IP22

Prevention of access to hazardous parts with fingers and protects the equipment inside the enclosure against ingress of solid foreign objects having a diameter of 12.5 mm and greater. Protected against vertically falling water drops when enclosure tilted up to 15°.

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All Rights Reserved. / Tous droits réservés. / Todos los derechos reservados. /
NuFACE® and Trinity® are registered trademarks of the Carol Cole Company DBA NuFACE.
Medical Device Safety Service GmbH