NŪFACE®

USER MANUAL

ENSURE YOUR BEST RESULTS

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Apply a mask-like layer of your Microcurrent Activator in sections as you lift. Push the ON/OFF button on your device and wait for the beep. Apply the Microcurrent Activator in the section you are about to lift and begin gliding the device until you hear the beep. Keep both spheres in contact with your skin during the treatment. Complete the same sequence in sections on your neck, cheeks, and forehead.

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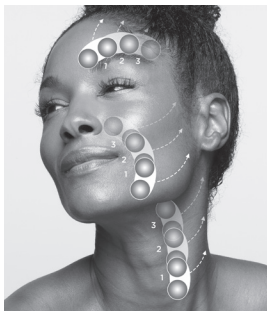
3) Finish / Terminer / Finaliza / Abschließen / Terminare / 마무리

Massage into skin any remaining Microcurrent Activator and finish by applying eye cream, moisturizer, and SPF.

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**2) Magnetic – Facial
Trainer Attachment**

Accessoire entraîneur facial magnétique
Accesorio magnético de entrenamiento facial
Magnetischer Aufsatz des Gesichtstrainers
Accessorio magnetico per l'allenamento facciale
마그네틱 - 페이스 트레이너 보조장치

3) Microcurrent Boost Button

Bouton Boost pour amplifier le microcourant
Elevador de la intensidad de la microcorriente
Taste zur Erhöhung des Mikrostroms
Pulsante Boost per potenziare la microcorrente
미세전류 부스터 버튼

4) Treatment Level Indicator

Indicateur de niveau de traitement
Indicador del nivel de intensidad
Anzeige der Behandlungsintensität
Indicatore del livello di trattamento
트리트먼트 수준 표시기

5) Dual Function Power Button: ON/OFF + Intensity Control Up/Down

Bouton d'alimentation à double fonction: March/Arrêt + Contrôle d'intensité haut/bas
Botón de encendido doble función: ENCENDIDO/APAGADO + Control de intensidad Arriba/Abajo
EIN-/AUS-Schalter und Intensitätsregelung auf/ab
Pulsante di alimentazione a doppia funzione: ON/OFF + Controllo d'intensità Sopra/Sotto
이중 기능 전원 버튼: 기기 켜기/끄기 + 강도 조정 감/약

6) Power Adapter (Not Shown)

Adaptateur (non illustré)
Adaptador de corriente (no se muestra)
Netzteil (nicht abgebildet)
Adattatore (non illustrato)
전원 어댑터(표시되지 않음)

NŪFACE[®] MINI+

GETTING TO KNOW YOUR DEVICE
SE FAMILIARISER AVEC L'APPAREIL
CONOCIENDO SU DISPOSITIVO
DAS GERÄT KENNENLERNEN
CONOSCERE IL DISPOSITIVO
기기에 대해 알아보기

1) Spheres (Electrodes)

Sphères (électrodes)
Esferas (electrodos)
Kugeln (Elektroden)
Sfere (elettrodi)
구체(전극)

2) Treatment Level Indicator

Indicateur de niveau de traitement
Indicador del nivel de intensidad
Anzeige der Behandlungsintensität
Indicatore del livello di trattamento
트리트먼트 수준 표시기

3) Dual Function Power Button: ON/OFF + Intensity Control Up/Down

Bouton d'alimentation à double fonction: March/Arrêt + Contrôle d'intensité haut/bas
Botón de encendido doble función: ENCENDIDO/APAGADO + Control de intensidad Arriba/Abajo
EIN-/AUS-Schalter und Intensitätsregelung auf/ab
Pulsante di alimentazione a doppia funzione: ON/OFF + Controllo d'intensità Sopra/Sotto
이중 기능 전원 버튼: 기기 켜기/끄기 + 강도 조정 강/약

4) Power Adapter (Not Shown)

Adaptateur (non illustré)
Adaptador de corriente (no se muestra)
Netzteil (nicht abgebildet)
Adattatore (non illustrato)
전원 어댑터(표시되지 않음)

The subject device produces low levels of microcurrent which is discharged through the two built-in chrome-plated spherical probes as they come in contact with the skin. Prior to using the subject device, the selected surface of the skin to be treated is first prepped with a NuFACE® Microcurrent Activator to moisturize the skin and ensure optimal microcurrent conductivity. The two spherical probes are then gently glided over the skin to deliver the electrical impulses per the selected treatment setting. Microcurrent is an aesthetic modality providing electric current in millionths of an ampere and has the ability to increase the appearance of facial contouring and the firmness of the skin and muscles.

INTENDED USE

The NuFACE Trinity+ and Mini+ devices are intended for facial and neck stimulation and are indicated for over-the-counter cosmetic use.

IMPORTANT INFORMATION

- Use the device only as described within this User Manual.
- The risks and dangers of using the device in any way that does not follow this User Manual are unknown.
- Register your new device at www.NuFACE.com.

CONTRAINDICATIONS FOR USE

- **DO NOT** use this device if you have a history of epilepsy or seizures.
- **DO NOT** use if you have had recent surgery or other surgical procedures.
- **DO NOT** use if you have any active cancer or over any suspicious or cancerous lesion.
- **DO NOT** use if you are pregnant. The safety of electrical stimulation during pregnancy has not been established.
- **DO NOT** use if you have an electronic implanted device such as a defibrillator, neurostimulator, pacemaker, or ECG Monitor.
- If you use a body-worn medical device, such as an insulin pump, consult your healthcare provider.
- **DO NOT** use on infants and children under the age of 18.
- Consult your healthcare provider if any of the contraindications above apply to you.

WARNINGS

- If you are under the care of a physician, consult your physician before using this device.
- This device is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Apply stimulation only to normal, intact, clean healthy skin. **DO NOT** use over swollen, infected, or inflamed areas or skin eruptions (e.g., phlebitis, thrombophlebitis, broken capillaries, varicose veins, etc).
- The product may **NOT** be safe for use on the neck in individuals with a known history of carotid artery disease, stroke, or transient ischemic attack (TIA), carotid stenosis, unstable blood pressure, easy fainting, thyroid disease, or prior surgical procedures in the neck area.
- **DO NOT** operate the device near electronic monitoring equipment (e.g., heart monitors, heart tracing or in close proximity (e.g. 1m or just over 3 feet) to shortwave, microwave, or electromagnetic therapy equipment as it may cause the medical equipment to malfunction. Such equipment may produce instability in device output and cause burns at the stimulation site.
- **DO NOT** disassemble the device as this may cause damage, malfunction, electrical shock, fire, or injury.
- Keep the device and related electrical accessories away from direct contact with water (including baths, showers, and sinks).
- **DO NOT** use across your chest as the introduction of electrical current into the chest may cause rhythmic disturbances to your heart, which could be lethal.
- Avoid applying excessive pressure when using the device on the anterior (front area) of the neck.
- **DO NOT** use over the midline of the neck. Care should be taken to minimize treatment directly over the midline of the neck, specifically avoiding the thyroid (see illustration).
- Keep the device out of the reach of children and pets.
- **DO NOT** drop or insert any object(s) into openings within the device.

Do not use
near the midline
of the neck

WARNINGS (CONTINUED)

Consult your healthcare provider before using this device:

- If you have any pre-existing or potential health conditions that may be impacted by use of the device.
- If you have painful areas.
- If you suspect or have a neuromuscular disorder.
- If you have a spinal cord injury.
- If you have any underlying skin conditions.
- If you have undiagnosed or known pain syndromes.
- If you have a tendency to bleed excessively following an injury.
- If you experience constant headaches and/or painful sensations.
- If you have suspected or diagnosed heart disease.

DO NOT apply stimulation or use the device:

- When there is a tendency to bleed excessively following acute injury or breaks in bones.
- On the eye area (circular muscle within the orbital rim).
- On the temples, the breast or groin area.
- Over painful areas of your skin or in areas that lack normal sensation.
- While bathing or showering.
- At extreme temperatures.
- If you suspect the device is damaged or not functioning properly.
- While driving or operating machinery.

Discontinue use of the device and contact your healthcare provider if you experience any of the following:

- Persistent pain or discomfort.
- Discomfort or any difficulty breathing or swallowing.
- Negative or harmful reaction(s) or injury.
- Symptoms of hypertension (dizziness or lightheadedness).
- Constant headaches and/or painful sensations.
- Skin irritation or burning.
- Unpleasant sensations on skin during use.
- Consistent flashing lights when without the use of your device.

If you experience any negative reactions while using the device and/or the microcurrent activator, discontinue use immediately and contact your healthcare provider.

PRECAUTIONS

- The NuFACE device is designed for cosmetic use only.
- The device and attachments should be operated, transported, and stored in the manner specified in this manual.
- Failure to use and maintain the device and corresponding attachments in accordance with the instructions outlined in this manual will void the warranty and may result in product damage or physical injury.
- The long-term effects of transcutaneous electrical stimulators are unknown.
- Apply stimulation only to normal, intact, clean, healthy skin.
- Avoid applying excessive pressure when using the device.
- Skin irritation or hypersensitivity may occur due to the electrical stimulation and / or use of the NuFACE® Microcurrent Activator.
- A microcurrent sensation may occur while using the device. Using a lower intensity and/or additional activator may be needed to reduce the sensation.
- Skin may become dry during treatment in which Microcurrent Activator may need to be reapplied.
- A perception of flashing lights may occur if device is used near or around the eye region.
- **DO NOT** use the device if you see any signs of damage or leaks from the device.
- **DO NOT** remove the battery of the device.
- **DO NOT** store or leave the device under the sun or on a hot surface.
- Follow all local regulations and waste management guidelines for safely disposing devices containing batteries.

PRECAUTIONS (CONTINUED)

- Only use the external power supply (marked NuFACE) provided with the device with the following specifications: 5VAC 1A. The use of accessories, transducers, and cables other than those sold by NuFACE may result in increased emissions, decreased immunity of the system, or injury.
- Always position the power adapter so it is accessible in case it needs to be disconnected from the wall outlet.
- Children should be supervised to ensure that they do not play with the device.

The NuFACE Trinity+ and Mini+ devices are intended for cosmetic use on the face and neck as described. Any consequence resulting from misuse, application of the device to other parts of the body or face that are not approved for use, connection to improper voltage sources, dirty spheres, expired or non-NuFACE conductive alternative, or any other improper application is neither the responsibility of NuFACE nor its affiliates.

The device has been designed to provide freedom from unacceptable risk when used under normal conditions and reasonably foreseeable hazards. In the event of complete loss of treatment function, the device does not pose a risk to the user.

TREATMENT AREAS AND TECHNIQUE

The instructions below are designed to provide users with the most effective results when followed correctly.

TREATMENT

- The device is designed to stimulate large surface areas of the face and neck, such as the jawline, cheekbones, and forehead. **DO NOT use the device directly over the midline of the neck.**

TECHNIQUE

- Using light to medium pressure, lift one side of the face and neck at a time using a gliding motion.

IMPORTANT: Always ensure that both spheres are in contact with the skin during treatment.

- For more advanced treatment options and holds, download the NuFACE App, or visit our website at www.NuFACE.com.

BEFORE USING YOUR NEW DEVICE

Follow these important steps before using your new device:

1) CHARGE YOUR DEVICE

Trinity+: Charge your Trinity+ device for ten (10) hours before the first use or until all three (3) LED indicator lights are lit solid.

- Remove your Trinity+ device and external power adapter from the box.
- Connect the external power adapter and plug the adapter into the wall outlet, keeping it away from any source of water.
- Place your Trinity+ device on the charging cradle.
- The lights will increasingly illuminate, indicating that your device is charging.
- Your device is fully charged when all the all the lights stay illuminated.

Mini+: Charge your Mini+ device for four (4) hours before the first use or until all three (3) LED indicator lights are lit solid.

- Remove your Mini+ device, and external power adapter from the box.
- Connect the external power adapter to the Mini+ device and plug the adapter into the wall outlet, keeping it away from any source of water.
- The lights will increasingly illuminate, indicating that your device is charging.
- Your device is fully charged when all the all the lights stay illuminated.

The device will NOT turn on when it is charging.

BEFORE USING YOUR NEW DEVICE (CONTINUED)

2) TEST YOUR SENSITIVITY

IMPORTANT: Test your skin to make sure you are not sensitive to the NuFACE Microcurrent Activator.

- Rub a small amount of Microcurrent Activator on your forearm and leave for **15-30 minutes** before washing it off.
- If an allergic reaction (rash, redness, etc.) develops, **DO NOT continue application of the Microcurrent Activator.** Contact your healthcare provider if the allergic reaction persists. Go to www.NuFACE.com to view other alternate NuFACE activators.

HOW TO USE YOUR DEVICE

Get ready to Tone, Lift, and Contour.

The instructions below are designed to provide the most effective results when followed correctly.

1) CLEANSE

- Cleanse your face and neck with an oil-free cleanser. Men should shave before treatment — facial hair can interfere with the conductivity of the spheres.
- **Note: Any cream, lotion, self-tanner, and/or oil-based products on the skin can block the microcurrent.**

2) ACTIVE

IMPORTANT: Always apply a NuFACE Microcurrent Activator when using your device.

- Apply a mask-like layer of the NuFACE Microcurrent Activator section-by-section as you treat.
- Reapply activator if skin becomes dry during treatment in order to continue safely using device.

3) PREPARE YOUR DEVICE

- **Trinity+:** Remove the power adapter from the device
- **Mini+:** Remove the power adapter from the device

4) SELECT YOUR TREATMENT MODE

- LEDS will indicate each treatment mode
- Press power button once, Instant Lift Mode ●○○○
- Press power button twice, Pro Toning Mode ○●○○
- Press power button three times, Skin Tightening Mode ○○○●

5) POWER ON YOUR DEVICE

- Press and hold the power button for 2 seconds
- An ascending audible beep is emitted, the energizing hum will begin, and the lights will illuminate indicating that your device is ON.
- Adjust the intensity level to your comfort by pressing the ON/OFF button on the front of your device.
- As you increase or decrease intensity, an audible beep is emitted, and the number of illuminated lights will change.

DEVICE INDICATORS

INTENSITY LEVEL	INDICATOR LIGHT	
Minimum Intensity	First Light Illuminated	●○○○
Medium Intensity	First & Second Light Illuminated	●●○○
Maximum Intensity	All (three) 3 Lights Illuminated	●●●●

HOW TO USE YOUR DEVICE (CONTINUED)

DEVICE AUDIBLE SOUND	DEVICE INDICATION
Ascending Beeps & Hums	Powered ON
Descending Beeps	Powered OFF
Single Beep	Intensity Adjusted
Dual-Tone Beep	End of Glide Interval (Tone occurs only during treatment)
Low-Tone Beep	Battery Needs To Be Charged (Tone occurs only during treatment)

6) BEGIN TREATMENT

- Place the device in direct contact with the treatment area.
- Always keep both spheres in contact with the skin during treatment to keep microcurrent flowing.**
- Using light to medium pressure, begin treatment.
- The device emits a dual tone beep every five (5) seconds to indicate when to move to the next treatment area.

PRO TIP FOR TRINITY+

During treatment, give stubborn areas an extra boost of NuFACE microcurrent energy by turning ON the **Boost button**.

- The energizing hum will increase, indicating that the boost is ON.

IMPORTANT: If you feel the electrical current and / or experience a tingling sensation on your skin, do the following:

- Add more NuFACE Microcurrent Activator to the area you are treating. The device will NOT contract the muscles when used as instructed.
- Lower the intensity on the device. Lowering the intensity may also reduce or eliminate this sensation.

7) FINISH TREATMENT

- Turn **OFF** the device after treatment.
- After treatment is complete and device is turned off, massage remaining NuFACE Microcurrent Activator into skin followed by eye cream, moisturizer, and SPF.

The device will automatically shut off after 20 minutes of use.

Go to the section '**CARING FOR YOUR DEVICE**' for instructions on cleaning and maintaining your device.

TREATMENT FREQUENCY

- Use 5 minutes a day.
- For the first 60 days, use 5 days per week.
- After 60 days, use 2-3 times per week for maintenance.
- Allow 24 hours between treatment sessions.

For more advanced treatment options and holds, download the NuFACE App, or visit our website at www.NuFACE.com.

RESULTS

Results may vary based on skin type and skin condition. Best results are achieved with regular use.

CARING FOR YOUR DEVICE

Cleaning:

- Slightly dampen a soft, lint-free cloth with water or rubbing alcohol to clean your device after every treatment.
- **DO NOT submerge your device in water at any time.**

Charging:

- Before first use, charge your **Trinity+ device for ten (10) hours or until all 3 LEDs are lit solid.**
- Before first use, charge your **Mini+ device for four (4) hours or until all 3 LEDs are lit solid.**
- The device will **NOT** turn on when it is charging.
- For continued optimal battery life, it is recommended to charge your device regularly.
- A fully charged device should last for approximately 2 weeks of use based on daily usage of 5 minutes.
- Keep the product out of the reach of children when charging.
- Keep the product away from high temperature, open flame and water sources when charging to prevent short circuiting or damaging the battery.
- Plug the charger into a wall outlet away from water, and within a location that can be accessed and disconnected easily.

Battery:

For ideal battery life, we recommend charging your device when the low battery indicator (1 LED light only) is displayed on your device and the indicator tone starts to descend.

Check Battery / Light Indicators:

To check your device battery level, **hold the power button down on your device for one (1) second.** The battery level will be illuminated through the three (3) LEDs:

LED LIGHT INDICATOR		BATTERY POWER
NO LED Lights	○○○	No Battery Power Trinity+: A full charge cycle of ten (10) hours is recommended before use. Mini+: A full charge cycle of four (4) hours is recommended before use.
1 LED Flashing	⦿○○	Critical Battery. Not enough charge to start a treatment.
1 LED	●○○	Low Battery
2 LEDs	●●○	Medium Battery
3 LEDs	●●●	Full Battery

Power Adapter:

- Only use the power adapter (marked NuFACE) provided with the device containing voltage 5VAC 1A.
- Connect the power adapter in a location that can be easily accessed and disconnected from the outlet.
- If your power adapter is damaged or missing, visit **www.NuFACE.com** to purchase a replacement.

Storage:

The device is intended to be stored and used at room temperature. Use caution prior to using the device if it has been stored or not used for an extended period of time. Inspect the device and do not use it if there are any signs of internal battery leaks.

Disposal:

IMPORTANT: If disposing the device, do not put into a fire. Dispose of the device and its batteries per local regulations.

TROUBLESHOOTING

Device Not Turning On:

1) Check to see if your device is fully charged.

- If the device does not turn ON, the battery may be drained.
- A full charge cycle of **ten (10) hours for Trinity+ and four (4) hours for Mini+** may be required for the device.
- For additional information go to the "Battery" section in this user guide.

2) For Trinity+ device only: Check to see if the attachment is properly seated and fits securely to the Trinity+ device.

IMPORTANT: The Trinity+ device will NOT turn on if the attachment is not correctly positioned and secured in place.

Resetting Your Device:

Trinity+:

- To reset your device back to factory settings, replace with plug your Trinity+ device to the power supply.
- Hold down the power button for **10 seconds**.
- The 1st and 3rd LEDs will blink three **(3) times**.
- Once the LEDs stop blinking you may start using your device.

Mini+:

- To reset your device, plug your Mini+ device to the power supply.
- Hold down the power button for **10 seconds**.
- The 1st and 3rd LEDs will blink three **(3) times**.
- Once the LEDs stop blinking you may start using your device.

DO NOT attempt to troubleshoot by disassembling the product as this will void the warranty and may cause injury.

CUSTOMER SUPPORT

If you are still experiencing issues with your device, contact NuFACE Client Services at info@nuface.com or call **888-360-9171** (US Only).

WARRANTY

NuFACE® warrants that all NuFACE® devices and attachments are free from manufacture defects in materials and workmanship for a period of two (2) years from the date of your original purchase. Our obligation under this Warranty is limited to replacement of the product which will not extend the period of your Warranty. This Warranty does not cover damage covered by abuse, improper care, application of harmful chemicals, product that may have been altered in any way, loss of the device and product purchased from an unauthorized seller.

NO OTHER WARRANTY EXPRESSED OR IMPLIED, WHETHER OF FITNESS FOR A PARTICULAR PURPOSE OR OF MERCHANTABILITY OR OF ANY OTHER KIND, SHALL EXIST IN RESPECT TO SUCH PRODUCT.

For Devices Purchased in the US: For more information about your warranty or to initiate a warranty claim and receive a return authorization number, contact our Client Services Team via our website at www.NuFACE.com or by calling **888-360-9171**. A receipt must be shown from an authorized seller as proof of purchase.

For Devices Purchased Outside of the US, International Returns, and International Warranties:

For all product purchased outside of the United States, warranty and returns are per your place of purchase—spa, retailer, or distributor. Please consult your receipt or place of purchase for their return and warranty policy.

CUSTOMER SUPPORT (CONTINUED)

TECHNICAL INFORMATION

The NuFACE® device complies with the requirements of IEC 60601-1-2:2014 (EMC Collateral Standard) including the E-field susceptibility requirements at 10 volts per meter, at frequencies from 80 MHz to 27 GHz. However, even at this level of device immunity, certain transmitting devices (cellular phones, two-way radios, cordless phones, paging transmitters, RFID devices, etc.) emit radio frequencies that could interrupt device operation if operated in a range too close to the device. Users should be aware of possible radio frequency interference if portable devices are operated in close proximity to the device. It is recommended that transmitting devices such as cellular phones, two-way radios, cordless phones, etc. not be located within 10 feet / 3m of the device during treatment. In addition, it is recommended that such transmitting devices not be located within 10 feet / 3m of the device during Bluetooth® wireless communication with the NuFACE App, other than the cellular phone being used to communicate with the device.

The NuFACE® device complies with the requirements of IEC 60601-1-2:2017 for Radiated and Conducted Emissions in Group 1 Class B.

The NuFACE® complies with the requirements of IEC 60601-1-2:2014 for Harmonic distortion in Harmonic classification A.

Mini+ Max Current: 335uA RMS
Trinity+ Max Current: 600uA RMS
Max Unloaded Voltage: 28VDC

FCC COMPLIANCE

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.
- Changes or modifications not expressly approved by NuFACE could void your authority to operate the equipment.

Industry Canada Compliance

This device contains license-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's license-exempt RSS(s). Operation is subject to the following two conditions: 1. This device may not cause interference. 2. This device must accept any interference, including interference that may cause undesired operation of the device.

CAN ICES-3 (B) / NMB-3 (B) This Class B digital apparatus complies with Canadian ICES-003.

Trinity+
REF
30790
10790-PRO

Mini+
REF
20400

Trinity+ FCC ID: 2AGNA-10760, IC: 25861-10760
Mini+ FCC ID: 2AGNA-20400, IC: 25861-20400

Conforms to EN IEC/UL: Device - 60529, 60335-1, 60335-2-23, 60601-1, 60601-1-11, 60601-1-2, 60601-1-6, 60601-2-10, 62133-2, 62366-1.
Power Supply - 60950-1, 60068-2-70

SYMBOLS

One or more of the following markings may appear on the device, User Manual or packaging.

European Conformity

Comply with European Restriction of Hazardous Substances Directive 2002/95/EC

European Representative

Product must be disposed of in accordance with the WEEE directive

Part Number

FCC-approved equipment authorization mark

Operator's Manual;
Operating Instructions

Class II double insulated electrical device

Refer to Instruction Manual/Booklet

Do not use this device in a bathtub, shower or water-filled reservoir

Manufacturer

Detachable Power Supply Part

Power On/Off

IP22

Trinity+/Mini+: Prevention of access to hazardous parts with fingers and protects the equipment inside the enclosure against ingress of solid foreign objects having a diameter of 12.5 mm and greater. Protected against vertically falling water drops when enclosure tilted up to 15°.

Type BF applied part, which is equipment that provides a particular degree of protection against electric shock. Particularly regarding allowable leakage current and of the protective earth connection.

RECYCLING

- Do not throw away the product with the normal household waste, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.
- This product contains a built-in rechargeable battery, which must be removed before you discard the product. **Due to the complexity of removing the battery, for your safety, we strongly recommend not to remove the battery by yourself. Please contact the after-sales service center (number indicated on the Chinese package label) and take your product to an official collection point to have a professional remove the rechargeable battery.**

Removing the rechargeable battery

Special Notes:

- Only remove the rechargeable battery when you are ready to discard the product. Make sure the battery is completely empty when you remove it.
- Take any necessary safety precautions when you use tools to open the product and when you dispose of the battery.
- Do not connect the product to the supply main after you have removed the battery.



Steps for removing the battery (Trinity+):

- Remove the attachment from the Trinity+ device.
- Cut along the gap reveal with a rotary saw to break the weld around the circumference of the Trinity+ device.


RECYCLING (CONTINUED)

3. Separate the front enclosure from the rear enclosure.
4. Cut the four battery retainer tabs with a rotary saw to remove the battery retainer.
5. Use snips to CAREFULLY cut through the two battery tabs.
6. Remove the battery for recycling.


(Please see the illustrations for each step below)

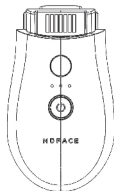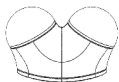
Step 1

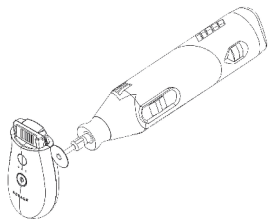Step 2

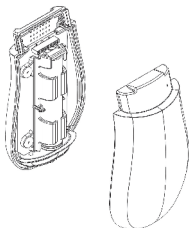Step 3

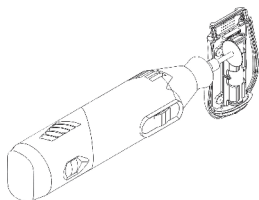
Step 4

RECYCLING (CONTINUED)

Step 5

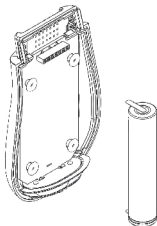
Step 6

Steps for removing the battery (Mini+):


1. Unplug the charging cable from the device.
2. Place a flathead screwdriver under the spheres and lift the spheres by twisting the screwdriver with moderate to high force.
3. Remove the spheres and seals from the device.
4. Use a rotary saw to cut the weld along the parting line around the circumference of the device.
5. Remove the front housing which contains the circuit board and battery.
6. Use the metal snips to CAREFULLY cut the battery tabs from the circuit board.
7. Remove the battery for recycling

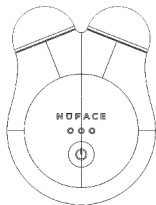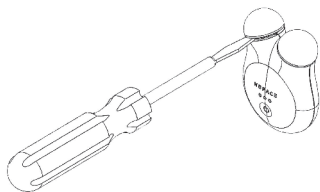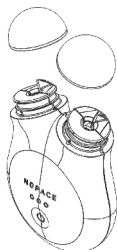
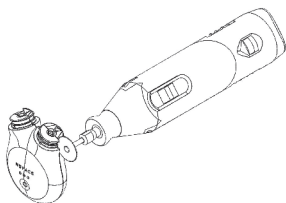
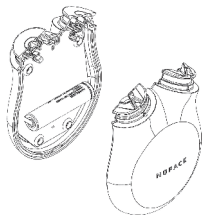
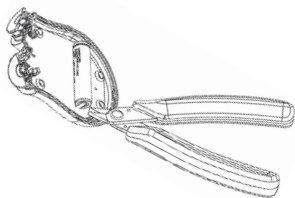
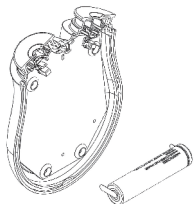
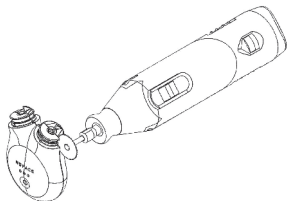
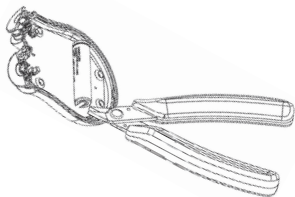
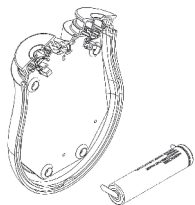
(Please see the illustrations for each step below)

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All Rights Reserved. / Tous droits réservés. / Todos los derechos reservados /
NuFACE® and Trinity® are registered trademarks of the Carol Cole Company DBA NuFACE.
Medical Device Safety Service GmbH