NÜFACE
TRINITY®

WRINKLE REDUCER
ATTACHMENT

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
NuFACE TRINITY® WRINKLE REDUCER ATTACHMENT

Accessoire de réduction des rides Trinity
Trinity Faltenentferner-Aufsatz
Trinity 红光嫩肤替换头 / Trinity 트리니티 링글 리두서

Note: Refer to NuFACE Trinity® User Manual for more device information.
Remarque : Référez-vous au manuel d'utilisation NuFACE Trinity® pour plus d'informations sur l'appareil.
Hinweis: Weitere Geräteinformationen finden Sie im NuFACE Trinity®-Benutzerhandbuch.
注意: 有关美容仪的更多信息请参阅《NuFACE Trinity® 使用说明书》, 참고: 자세한 장치 정보는 NuFACE Trinity® 사용 설명서를 참조하십시오.

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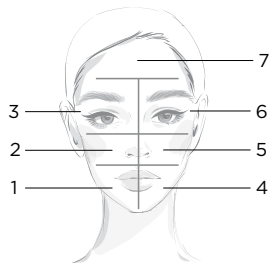1) NuFACE Trinity® Wrinkle Reducer Attachment / Accessoire de réduction des rides NuFACE Trinity®
NuFACE Trinity Faltenentferner-Aufsatz / NuFACE Trinity 红光嫩肤替换头 / NuFACE Trinity 트리니티 링글 리두서

2) Protective Shield and Skin Sensor / Surface de protection et capteur cutané
Schutzelement und Hautsensor / 防护罩/皮肤传感器 / 보호막/피부 센서

3) Power On/Off / Bouton Marche/Arrêt / Ein-/Ausschalttaste / 开/关按钮 / 전원 온/오프 버튼

4) Attachment Release Button (on back, not shown) / Bouton de déverrouillage de l'accessoire
(au dos, non illustré) / Aufsatz-Löseknopf (auf der Rückseite, nicht abgebildet)
替换头装卸按钮 (位于背面, 未显示) / 부가 장치 분리 버튼 (뒷면, 보이지 않음)

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SMOOTH - Place the NuFACE Trinity with Wrinkle Reducer Attachment directly on skin of one of the 7 targeted areas. Hold in place until you hear the beep (3 minutes), then move to the next area. Do not treat upper eyelid.

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HELLO FROM NuFACE®!

Congratulations on giving wrinkles the red light with your new Trinity Wrinkle Reducer Attachment.

We're excited for you to discover the incredible power of infrared and red light therapy. Simply switch out your Trinity Facial Trainer Attachment for your new Trinity Wrinkle Reducer Attachment and sit back as its LED lights help reduce the appearance of full-face wrinkles.

As always, don't forget to take before and after photos and tag us at @MyNuFACE. We'd love to follow your skincare journey!

Have an uplifting day,
Your NuFACE Family

INTENDED USE

The NuFACE Trinity Wrinkle Reducer is an over-the-counter hand-held device intended for the treatment of full-face wrinkles.

IMPORTANT INFORMATION

- The risks and dangers of using the Trinity Wrinkle Reducer in any way that does not follow this User Manual may be unknown. Use the Trinity Wrinkle Reducer only as described within this User Manual. The Trinity Wrinkle Reducer was not tested for use over the eye socket or eyelid so the risks are unknown.
- Your Trinity Wrinkle Reducer has a unique serial number located on the bottom of the attachment. Record your serial number and register your device online at www.MyNuFACE.com. A registration link can be found at the bottom of the webpage.
- In addition to guidance within this manual, also refer to your NuFACE Trinity Device manual.

CONTRAINDICATIONS FOR USE

- Do not use if you have a history of seizures, especially if you have photosensitive epilepsy.
 - Do not use on the neck if diagnosed with or suspect a thyroid problem.
 - Do not use if your skin is sensitive to light. People can be sensitive to light and not be aware of it. Certain drugs or cosmetics can cause skin sensitivity.
 - Do not use over any suspicious or cancerous lesion.
 - Do not use over topical or systemic steroids.
 - Do not use on infants and children under the age of 14.
- *Consult your healthcare provider if any of the contraindications above apply to you.

⚠️ WARNINGS

- If you are under the care of a physician, consult your physician before using this device.
- Discontinue using the device and consult your healthcare provider if you experience an adverse reaction from using the device.
- Apply only to normal, intact, clean, healthy skin.
- Consult your healthcare provider before use if you have any underlying skin conditions.
- Do not use flammable aerosol products or oxygen near the device.
- Do not operate Trinity Wrinkle Reducer in the presence of electronic monitoring equipment (e.g., heart monitors, heart tracing or in close proximity (e.g. 1 m) to shortwave, microwave or electromagnetic therapy equipment as it may cause the medical equipment to malfunction. Such equipment may produce instability in the Trinity Wrinkle Reducer.
- If the device begins to overheat, turn the device OFF immediately.
- Never operate the device if it's not functioning normally, or not operating properly due to damage caused by dropping, submersion in water or accidental cracks.
- Never drop or insert any object into any openings on the device.
- Device should not be used at extreme temperatures.
- Do not look directly into the Trinity Wrinkle Reducer's LEDs or shine in the eyes as this may cause harm.
- Overuse of this device may result in dryness of the skin.
- Keep all electrical appliances (including the Trinity Wrinkle Reducer) away from water (including baths, showers, and sinks).
- Do not use while driving, operating machinery, or during any activity which may put the user or the public at unnecessary risk of injury.
- Discontinue using the device immediately if you suspect the device to be malfunctioning, or not operating properly. Return the device to NuFACE for examination and repair.

⚠ WARNINGS (CONTINUED)

- Use of controls or adjustments, or the performance of procedures other than those specified herein, may result in injury.
- Do not disassemble the attachment as this may cause damage, malfunction, electrical shock, fire or injury. There are no user-serviceable parts inside. **Note: Disassembly voids the warranty.**

PRECAUTIONS

The Trinity Wrinkle Reducer is designed for cosmetic use only and for individuals in good health. If this is not the case for you, please consult your healthcare provider before use.

- If using the Trinity Wrinkle Reducer causes persistent pain or discomfort, discontinue use immediately.
- The Trinity Wrinkle Reducer is intended to be stored and used at room temperature. Do not store Trinity Wrinkle Reducer in the sun or on a hot surface. High temperatures may be damaging. The device should be operated, transported and stored in the manner specified in this manual. Failure to use and maintain the device in accordance with the instructions outlined in this manual will void the warranty.

The NuFACE Trinity Wrinkle Reducer is intended for cosmetic use on the face as described. Any consequence resulting from misuse, application to other parts of the body or face, connection to improper voltage sources, dirty surface, or any other improper applications is neither the responsibility of NuFACE nor its affiliates.

GETTING STARTED

FOLLOW THESE IMPORTANT STEPS BEFORE USING YOUR NEW ATTACHMENT.

- 1) Note:** The Trinity Wrinkle Reducer's unique technology utilizes flashing LEDs to address facial wrinkles. Begin treatment in a well-lit room if you find flashing lights bothersome.
- 2)** Firmly press the release button to remove the current attachment. When fastening a new treatment attachment, press the release button, guide the attachment in place, and then release the button to secure the attachment to your Trinity Device.
- 3)** The +/- keys are not active when the Trinity Wrinkle Reducer is attached. There is only one setting for this attachment.
- 4)** Now, test your skin to make sure you are not sensitive to light. To do so, place the Trinity Wrinkle Reducer on your forearm for 3 minutes and then check the area. If your skin turns red and the redness lasts for more than 2 hours, you may be sensitive to light therapy and should consult with your healthcare provider before continuing to use the Trinity Wrinkle Reducer.
- 5)** Before you start, remember to take a photo of your cleansed face in a well-lit room to document your results. After your first treatment, take an after photo and repeat the photos in the same well-lit room every 4 weeks to track your results.

DEVICE SOUND — INDICATES

Ascending beeps	Powered ON
Descending beeps	Powered OFF
Dual-tone beep	End of interval
Three beeps	Skin not in contact

TREATMENT AREAS

The Trinity Wrinkle Reducer is designed to address large surface areas of the face. The following instructions are designed to provide users with the most effective results when followed correctly.

HOW TO USE YOUR TRINITY WRINKLE REDUCER

GET READY TO GIVE WRINKLES THE RED LIGHT.

- 1)** Cleanse your face with an oil-free cleanser.
- 2)** Remove your device from the charging cradle and turn it on by pressing the ON/OFF button. An ascending audible beep is emitted notifying you that your device is ON, and the red LEDs will flash rapidly 3 times. Do not touch the Trinity Wrinkle Reducer shield until the unit has finished flashing.
- 3)** Divide the face into 7 sections (see diagram in Quick Start) and treat each section separately.
- 4)** Place the Trinity Wrinkle Reducer in direct contact with your clean, dry skin, over the area to be treated. Your Trinity Wrinkle Reducer has a skin sensor that will illuminate the red LEDs when in complete contact with your skin. Press firmly enough to keep the Trinity Wrinkle Reducer in direct contact with your skin. The device should be used on lines and wrinkles within the section, not specifically where the number appears on the photo. Treat each section for 3 minutes until the beep sounds.
- 5)** If the Trinity Wrinkle Reducer is not touching the skin, your device will emit 3 beeps notifying you to hold directly on your skin.
- 6)** The Trinity Device will automatically shut OFF after 21 minutes indicating you are now finished with your treatment.

Note: Do not immerse device in water at any time.

TREATMENT FREQUENCY

- Use no more than 3 minutes per area, up to 21 minutes a day for the entire face. Do not overlap areas.
- For the first 60 days, use 5 days per week.
- After 60 days, use 2-3 days per week for maintenance.
- Allow 24 hours between treatments.

RESULTS

Results may vary based on skin type and skin condition. It takes about 8 weeks to see the start of the reduction in the appearance of fine lines and wrinkles. Some users see results in a shorter period of time, while others take longer.

CARING FOR YOUR TRINITY WRINKLE REDUCER

CLEANING

Slightly dampen a soft, lint-free cloth with water or rubbing alcohol to clean your device and attachment. **Do not immerse device in water at any time.**

CARING FOR YOUR TRINITY WRINKLE REDUCER (CONTINUED)

TEMPERATURE SENSOR

It's normal for the Trinity Wrinkle Reducer to get warm during treatment. If the attachment's internal temperature exceeds 129.5°F, a temperature sensor will automatically turn the device off. Contact NuFACE client services at 888-360-9171 if this occurs.

TROUBLESHOOTING

If your Trinity Device is not turning on, please make sure all of the following are true:

- The Trinity Wrinkle Reducer is completely attached to the device and fits snugly. Make sure to firmly press the release button when attaching and detaching the Trinity Wrinkle Reducer.
- The Trinity Device battery is fully charged. If not, please charge your device for a full 12 hours.
- The Trinity Wrinkle Reducer is pressed completely to the skin, activating (or turning on) the LEDs. Remember, the Trinity Wrinkle Reducer has a skin sensor that only turns on the LED lights when in direct contact with the skin.

If your device is still not working, please call NuFACE Client Services at 888-360-9171. We're here to help!

WARRANTY

For Devices Purchased in the US: NuFACE warrants that the NuFACE® devices and attachments are free of Manufacture defects in materials and workmanship for (1) year from the date of purchase for the NuFACE Trinity® Device, (2) years for the Trinity PRO Device. Abnormal wear and tear, submersions and drops are not covered by this warranty. This warranty is void if the product has been subject to accident, misuse, abuse, improper maintenance or repair, unauthorized modification, or any other use not found in the printed instructions. **Must show receipt as Proof of Purchase.**

For more information on your warranty, visit
www.MyNuFACE.com/pages/warranty-returns

Procedure for Warranty Service: If warranty service is necessary you can download a Warranty Form on our website at: www.MyNuFACE.com/warranty or call Client Services at 888-360-9171 to receive further instructions. Please be sure to return all parts (device, cradle and charger) and include a completed Warranty Form.

NuFACE - Returns / Repairs
1325 Sycamore Ave, Suite A
Vista, CA 92081
USA
1-888-360-9171

For Devices Purchased Outside of the US, International Returns and International Warranties:

All product returns outside of the United States are per your place of purchase - spa, retailer or distributor. Please consult your receipt or place of purchase for their return and warranty policy.

Conforms to EN IEC/UL: 60601-1, 60601-1-6, 60601-1-11, 60601-2-10, 60601-1-2, 60601-1-4, 60335-1, 60335-2

SYMBOLS

One or more of the following markings may appear on the Trinity Device, User Manual or packaging.

European Conformity

European Representative

Part Number

Refer to Instruction Manual

Manufacturer

Type BF applied part, which is equipment that provides a particular degree of protection against electric shock. Particularly regarding allowable leakage current and of the protective earth connection.

Comply with European Restriction of Hazardous Substances Directive 2002/95/EC

Product must be disposed of in accordance with the WEEE directive

FCC-approved equipment authorization mark

Class II double insulated electrical device

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All Rights Reserved. / Tous droits réservés. / Alle Rechte vorbehalten. / 版权所有 / 모든 권리 보유.
NuFACE® is a registered trademark of the Carol Cole Company.
Medical Device Safety Service GmbH